

POST-OP KNEE BRACING



PREMIUM SIZED KNEE IMMOBILIZER

- Tri-laminate construction of canvas, foam and brushed nylon.
- Elastic patellar strap.
- Tapered Velcro straps with loop lock closures
- Sturdy posterior, medial and lateral stays

PART#	SIZE	THIGH CIRC.*	LENGTH
200023	Small	36-41cm	20"
200025	Medium	41-46cm	20"
200027	Large	46-51cm	20"
200028	XLarge	56-56cm	20"
200029	XXLarge	56-61cm	20"

* Measure thigh circumference 10cm above the knee




Cool Version

POST-OP KNEE BRACING – INNOVATOR® - COOL VERSION

The simple, versatile post-operative knee bracing solution

- Locks quickly at 0°, 10°, 20°, 30°, 45°, 60°, 75° & 90°

PRODUCT HIGHLIGHTS

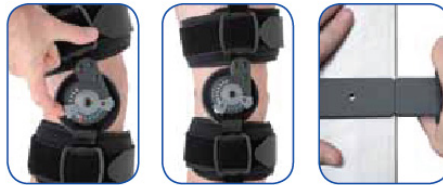

1. Simply press, rotate and release flexion/extension settings.
2. Low-profile hinge, with no exposed metal.
3. Break-away struts for easy, tool-free length adjustment.

PART#	FOAM	SIZE	THIGH CIRC.*	LENGTH
403100	Cool	Small	Up to 69cm	24"

Indication: - Post-op/post-injury ROM control or immobilization



POST-OP KNEE BRACING – EPO LITE

An economical solution for post-operative indications, the EPO Lite features the Quick-Set hinge and has reduced foam coverage for greater breathability and comfort. Universally Sized.

Features

- Quick-Set hinge for range of motion control.
- Reduced foam for breathability and comfort.

PART#

09360	EPO Lite, Short
-------	-----------------

*Circumference taken 6" (15 cm) above mid patella

Indication: - ACL, PCL, MCL, LCL injuries

- Tibial plateau fractures
- Osteochondral repairs
- Meniscal repairs
- Patella tendon repairs
- Condylar fractures
- HTOs (High Tibial Osteotomy)
- Acute sprains/strains of the knee

POST-OP KNEE BRACING

POST-OP KNEE BRACING - T SCOPE® PREMIER

Features - Performance

- Easy to use ROM hinge with quick-adjusting flexion and extension stops to control and limit knee flexion and extension.
- Telescoping calf and thigh sleeves for sizing a wide range of patients from 5' to 6'4" tall.
- Easy to use extension drop lock allows the brace to be locked out in 5 positions (in -10°, 0°, 10°, 20°, and 30°) of knee extension with the push of a button.
- Stiff telescoping bars and overall brace design limits knee range of motion.
- Accommodates both right and left leg procedures, reducing the need to stock additional products.

Easy to use

- Extremely quick application time and simple to adjust to enhance patient compliance and comfort.
- Easily and quickly customized for fit across a wide range of patient sizes.
- Trimmable straps for size customization.
- Strap lock clips keep straps positioned on bars and the hinges in the correct placement around the knee.
- Quick clip buckles allow for easy application and removal of the brace.

Comfort

- Comfortable padding throughout the brace, designed with patient compliance in mind.
- Brace is lightweight, weighing less than (1 kg) (standard version).
- BridgeTech Incision pad relieves pressure and provides added comfort around the incision site (optional accessory).
- Excellent overall suspension due to the design of the pads, strapping and overall brace.
- More flexible cuffs allow for better comfort and suspension of the brace.
- Near universal sizing, with the hinge bars extending from 17" to 27" to accommodate a wide range of patient anatomies.
- Standard version has reduced foam vs. other post-op braces for greater breathability.

PART#

08814	T Scope Premier 17"-27", fits up to a 30.5" thigh, universal Right or Left
73621	T Scope Premier Incision Pad

Indication: - ACL, PCL, MCL, LCL repairs/injuries

- Tibial plateau fractures
- Osteochondral repairs
- Meniscal repairs
- Patella tendon repairs
- Condylar fractures
- Sprains/strains of the knee
- HTOs (High Tibial Osteotomy)



POST-OP KNEE BRACING – T-SCOPE

The most versatile of all Post-Op braces, the T Scope features telescoping hinge bars that are adjustable from 17" to 27" along with friction locks to easily secure the brace in position once the desired length is selected. The T Scope also has quick clip buckles and strap lock clips that allow for easy fitting. Universally sized.

Features

- Telescoping hinge bars extendable from 17" – 27" (43.2 cm - 68.6 cm) accommodate a larger range of patient leg lengths.
- Lightweight - Weighs less than 37 oz. (1.05 kg).
- Universally sized (fits up to 30.5" thigh circumference).
- Strap lock clips keep straps positioned on bars.
- Quick clip buckle allows for easy on/off.
- Quick lock offers ability to lock out brace from -10° to 30° (with increments of 10°).
- Reduced foam for greater breathability.



PART#	DESCRIPTION	CE Marked
07714	T-SCOPE Post-Op	